

最薄纸张不覆膜

20X11cm

折

折

折

隔5CM折叠

50.00 mm

### WASH CARE

Only eye mask cover is washable. Do not allow the inner heating pad to get wet.

### WARNING

DANGER-To reduce the risk of burns, electric shock and fire, this product must be used in accordance with the following instructions.

- Read all instructions.
- Do not use while sleeping.
- Never pull this eye mask by the supply cord. Burns will result from improper use.
- Stop heating or switch to a lower temperature setting when reaching the desired temperature. Even comparatively low temperature can cause moderate temperature burn, check your skin frequently.
- The product is not to be used on or by an invalid, sleeping or unconscious person, or a person with poor blood circulation unless carefully attended. The elderly and children should use this product under supervision.
- People who have cutaneous sensory disorder should not use this product.
- Save these instructions.

## 尾页



### Registration for Warranty

Your new Aroma Season product comes with one year long warranty against manufacturing defects. Register the 12 months quality warranty for your order via this link: [www.aromaseasoninc.com](http://www.aromaseasoninc.com). We will verify and send you a confirmation message.



### Contact Us

For support, warranty, replacements or returns, please contact us at [info@aromaseasoninc.com](mailto:info@aromaseasoninc.com), including your order number and website where you placed order, then we will do whatever we can to make things right.

### REGISTRATION FOR WARRANTY



[www.aromaseasoninc.com](http://www.aromaseasoninc.com)

## 首页

## AROMA SEASON® ELECTRIC HEATED EYE MASK

AROMA SEASON INC  
Email: [info@aromaseasoninc.com](mailto:info@aromaseasoninc.com)  
Address: 1150 S MILLIKEN AVE. #1061  
ONTARIO CA 91761  
MADE IN CHINA

正面


### Owner's Manual

Thank you for purchasing Aroma Season Electric Heated Eye Mask.  
For optimum operation, please read this manual carefully before use.

### SPECIFICATIONS

Item Name: Electric heated eye mask	Material: Polyester
Model No.: AM-HE820-PS	Dimension: 8.5*3.5 inches
Input: DC 5V/1A-2A	Output: DC 3-5V/Max 1A
Temperature Mode: High/Medium/Low	
Present Timer: 20mins/40mins/60mins	

### HOW TO USE

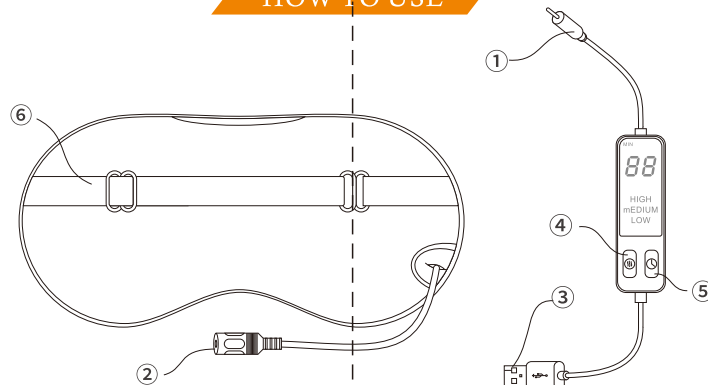

- ① DC plug on controller
- ② DC socket on eye mask
- ③ USB plug
- ④ Temperature mode button
- ⑤ Timing button
- ⑥ Adjustable strap

1. Insert the DC plug① on controller into the DC socket② on eye mask.
2. Insert the USB plug③ on controller into power source.
3. Press the <⏪> button④ on controller to start heating.
4. Press the <⏩> button④ to select a temperature mode(High/Medium/Low), and press the <⌚> button⑤ to select the heating time.
5. Spray 2-3 drops of warm water on eyes side of the eye mask for moist heating if needed, and adjust the strap⑥ to wear the eye mask around your head comfortably.
6. Long press the <⏪> button④ for 2 seconds to turn it off. The eye mask will stop working automatically when timing is over.

(Note: There is no battery inside the eye mask. Keep it power on when using it.)

反面